



World Health
Organization



Module III
Packaging of Infectious Substances

August 2011



World Health
Organization

Overview

- ☣ Triple Packaging
- ☣ Packaging Category A Infectious Substances (PI620)
- ☣ Packaging Category B Biological Substances (PI650)
- ☣ Packaging Exempt Human/Animal Specimens
- ☣ Packaging Limitations
- ☣ Packaging Exercise



World Health
Organization

How are Category A, Category B, and
Exempt Human/Animal Specimens
packaged?




World Health
Organization

Triple Packaging

Triple packaging provides three layers of containment to protect the substances being shipped. These layers are primary, secondary, and outer containers.

The following diagram shows the basic concept of triple packages.

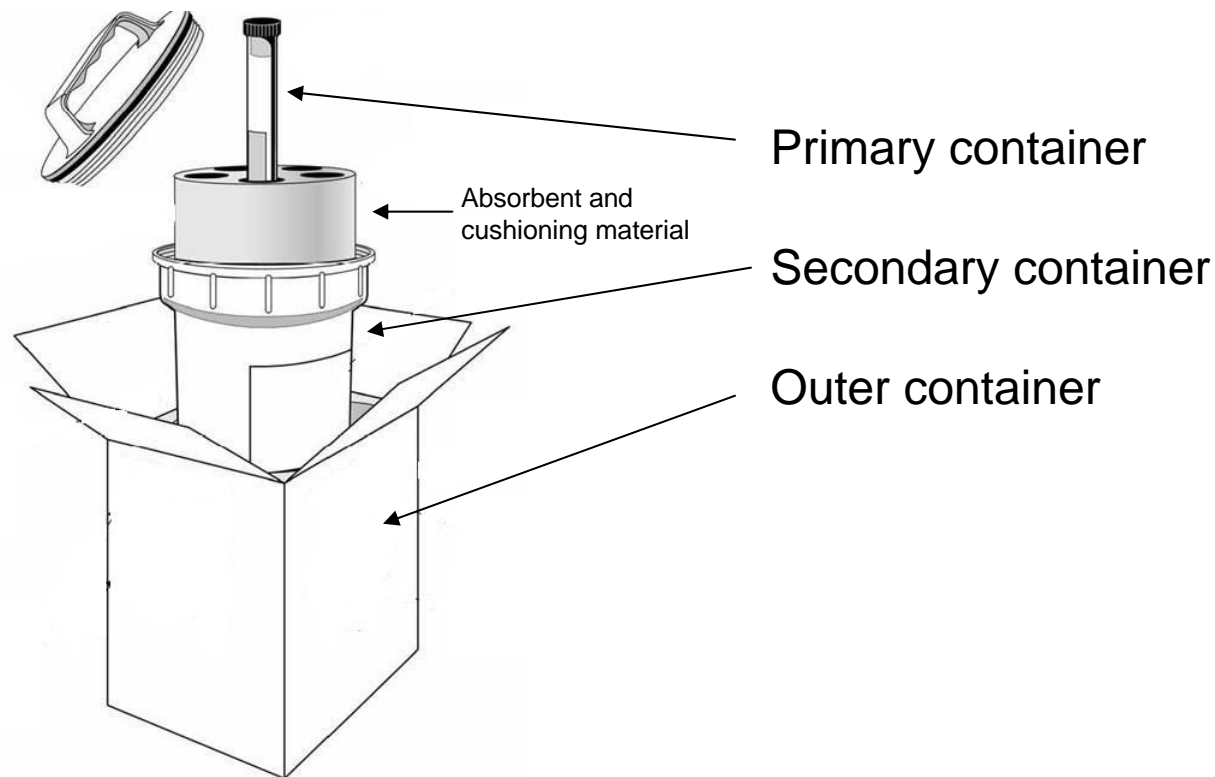

Primary container

Absorbent and
cushioning material

Secondary container

Outer container



**World Health
Organization**

What are the requirements for
packaging Category A substances?

Category A Packaging Requirements


- ☣ Primary container is leakproof
- ☣ Secondary container is leakproof
- ☣ Outer container is rigid
- ☣ UN specification marking:
 - ☣ Pressure tested at 95 kPa
 - ☣ Drop tested from 9 m
 - ☣ Puncture tested at 7 kg
 - ☣ Stacking tested
- ☣ Shipper must be trained





World Health Organization

Category A Packaging Requirements

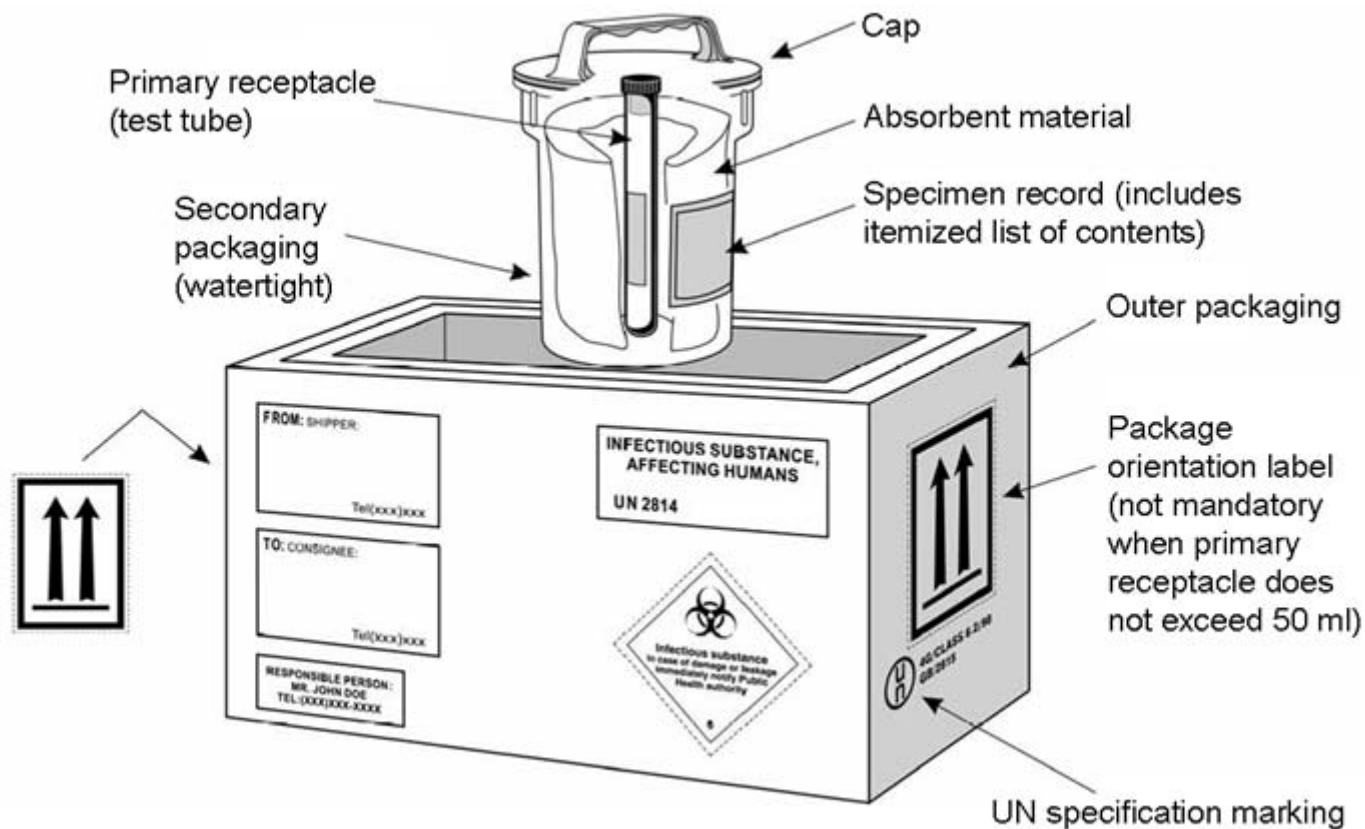



**World Health
Organization**

What are the requirements for
packaging Category B substances?

Category B Packaging Requirements

- ☣ Primary container is leakproof
- ☣ Secondary container is leakproof
- ☣ Either the primary or the secondary container must be pressure tested at 95 kPa
- ☣ Either secondary or outer container is rigid
If the shipment is transported by air, the outer container must be rigid.
- ☣ Drop tested from 1.2 m




World Health
Organization






World Health
Organization

Category B Packaging Requirements

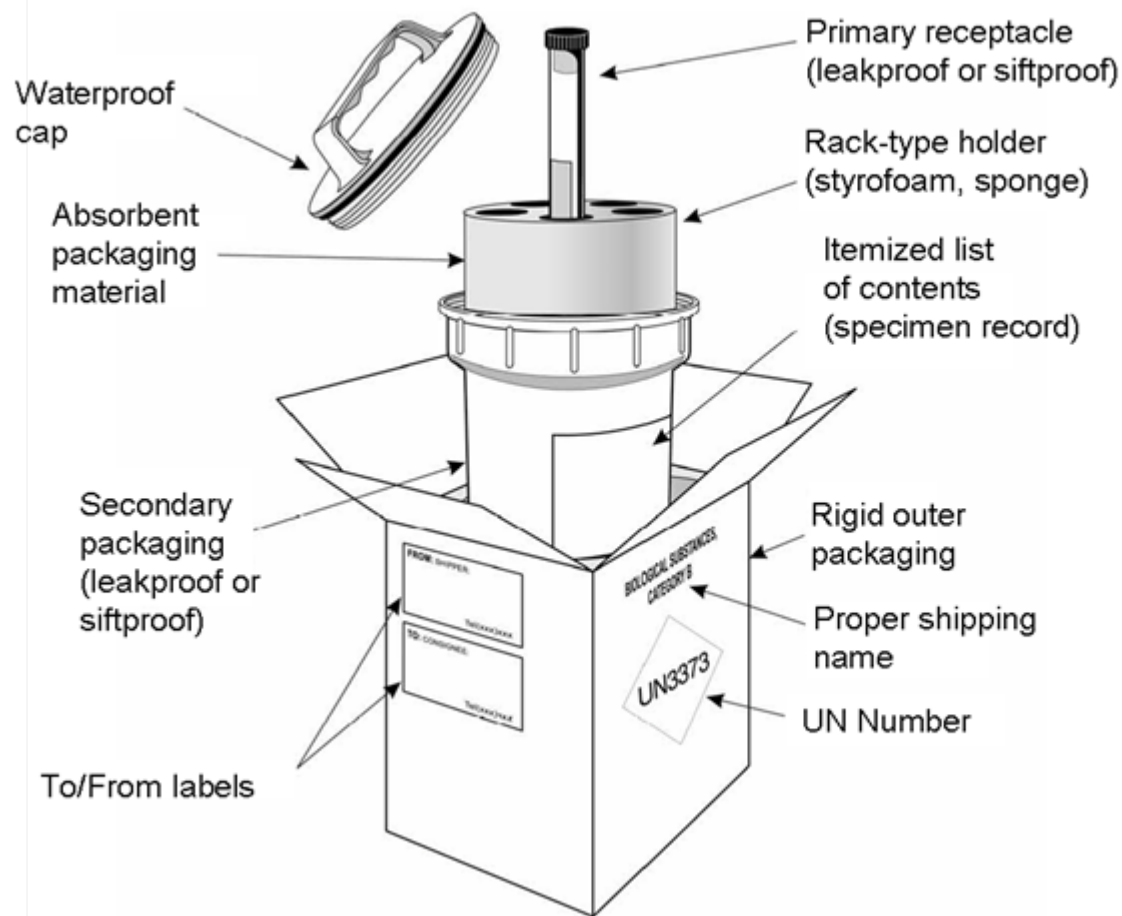

Category B Packaging Requirements

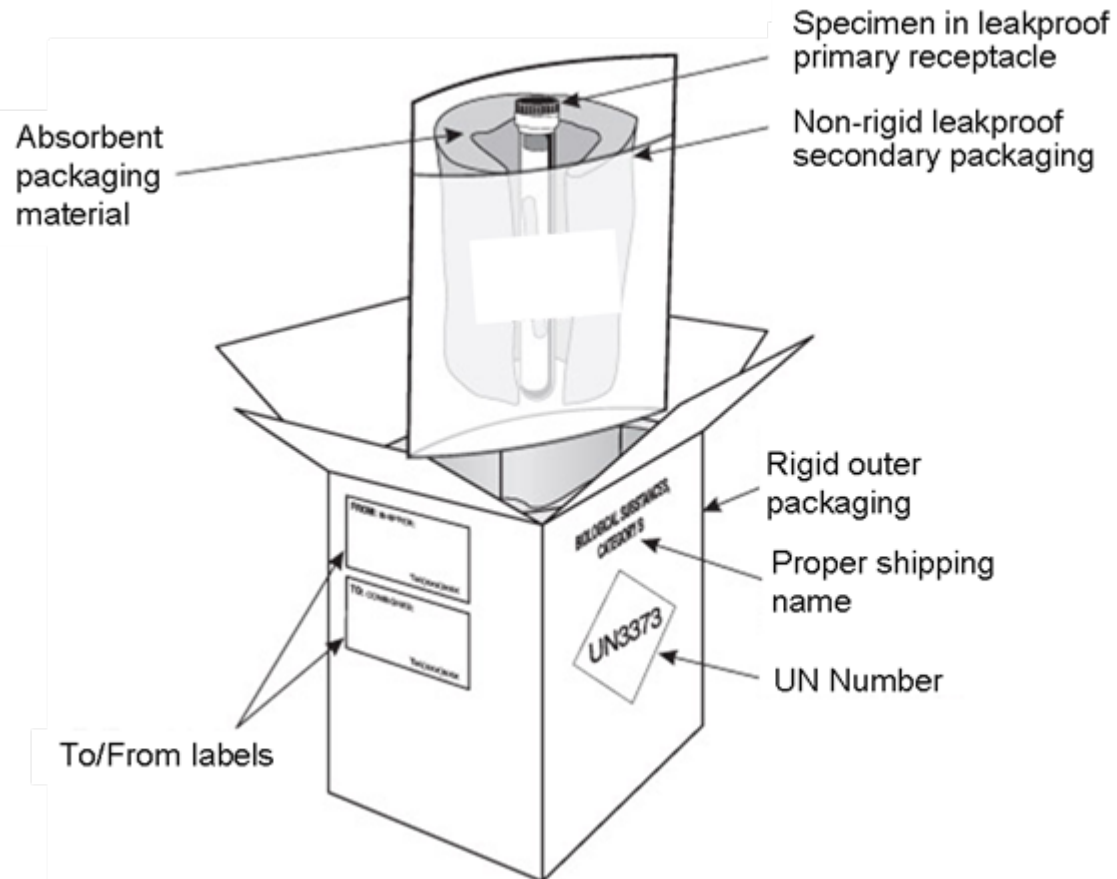



**World Health
Organization**

What are the requirements for
packaging Exempt Human/Animal
Specimens?

Packaging Requirements for Exempt Human/Animal Specimens


- ☣ Primary container is leakproof
- ☣ Secondary container is leakproof
- ☣ Outer container must be of adequate strength



World Health
Organization



Exempt Packaging Requirements

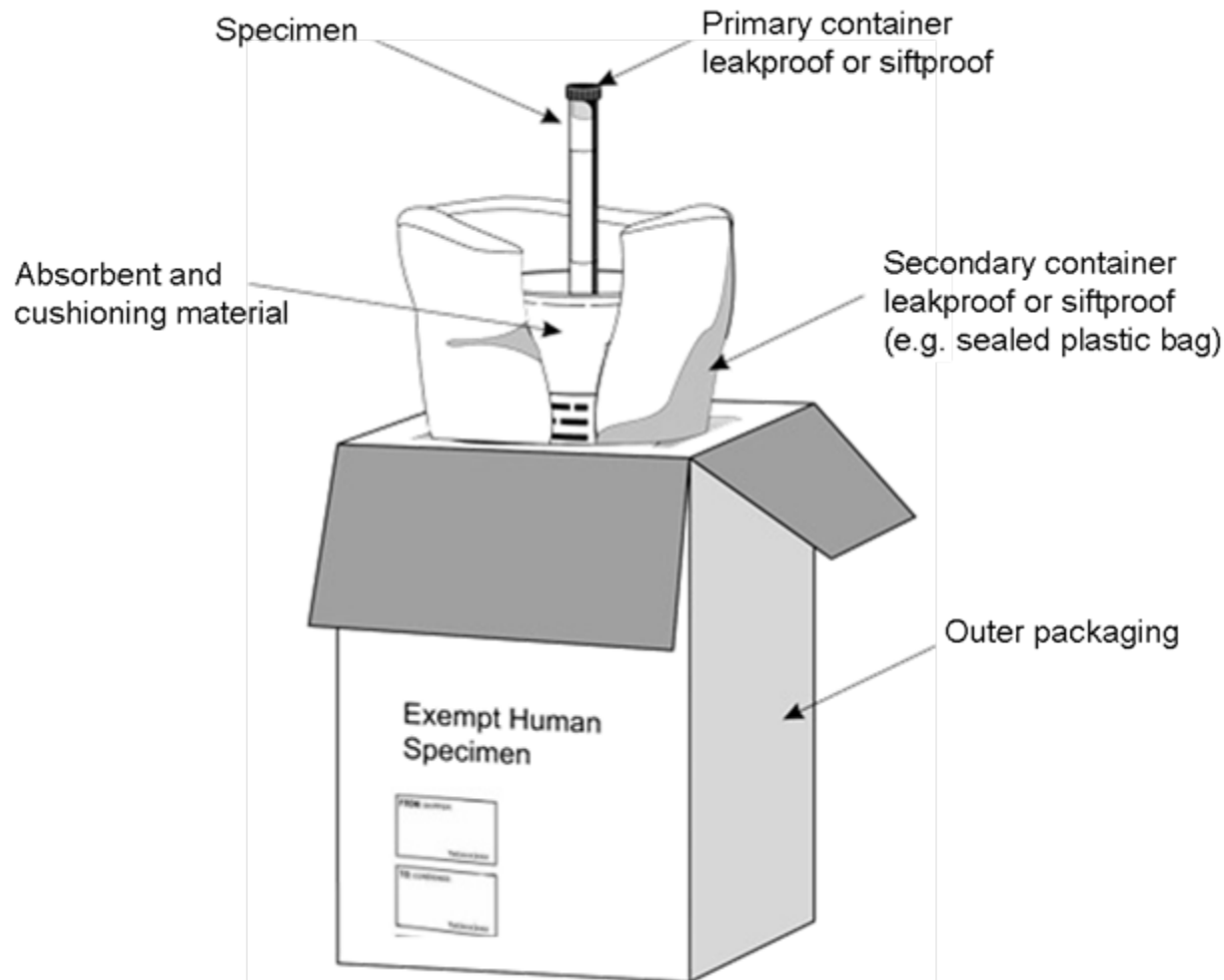



**World Health
Organization**

What are the limitations for shipping
infectious substances by air?

Limitations for Shipping

Weight and volume (excluding overpacks) Module II slide 16 columns 13 & 15

Category A Substances

Maximum 50 ml or 50 g per package for passenger aircraft

Maximum 4 litres or 4 kg per package for cargo aircraft

Category B Substances

Maximum 4 litres or 4 kg per package for passenger or cargo aircraft

Maximum 1 litre per primary container for passenger or cargo aircraft

Multiple samples

Multiple primary containers can be placed in the same package and must be individually wrapped or separated to prevent contact

Package dimensions

PI620 packages:

The smallest external dimension shall not be less than 10 cm

PI650 packages:

At least one surface of the outer packaging must have a minimal dimension of 10 cm x 10 cm





**World Health
Organization**

Can packages shipped to the same
destination be consolidated?



Overpacks

- ❖ Overpacks are enclosures used by a single shipper to contain one or more packages shipped to a single destination and to form one unit for convenience of handling, stowage or refrigeration during transport.
- ❖ Examples of overpacks are a number of packages either:
 - ❖ Placed or stacked on to a load board such as a pallet and secured by strapping, shrink wrapping, stretch wrapping, or other suitable means; or
 - ❖ Placed in a protective outer packaging such as a box or crate; or
 - ❖ Placed into an insulated chest and surrounded with refrigerants such as ice, dry ice or gel pads



Overpacks

- ☣ Overpacks are useful to ship samples that exceed the individual package quantity limits

For example:

- ☣ A shipper needs to send 100 ml of a Category A substance on a passenger aircraft. However, the quantity limit is 50 ml per package
- ☣ The sample could be split, divided into two packages of 50 ml, and then consolidated into a single overpack
- ☣ There are no specifications or testing requirements for overpacks
- ☣ However, any package placed in an overpack must comply with the regulations

Overpacks





World Health
Organization

Exercise III

Packaging Infectious Substances

Provided Supplies

- ☣ Gloves
- ☣ 2 Category A samples in primary containers
- ☣ Secondary container
- ☣ Rubber bands
- ☣ Absorbent material
- ☣ Inner lining
- ☣ Cushioning material
- ☣ Outer container
- ☣ Disinfectant
- ☣ Itemized list of contents



World Health
Organization



World Health
Organization

SOP for Category A Package

Step	Practice
1	
2	
3	
4	
5	
6	
7	
8	
9	
10	
11	
12	
13	



**World Health
Organization**

Questions and comments