



World Health  
Organization

## DISSOLUTION TEST FOR SOLID ORAL DOSAGE FORMS

### Draft proposal for revision in *The International Pharmacopoeia*

(June 2020)

*DRAFT FOR COMMENTS*

Please send any comments you may have on this draft working document to **Dr Herbert Schmidt**, Technical Officer, Norms and Standards for Pharmaceuticals, Technical Standards and Specifications ([schmidth@who.int](mailto:schmidth@who.int)) by **31 August 2020**.

Working documents are sent out electronically and they will also be placed on the WHO Medicines website ([http://www.who.int/medicines/areas/quality\\_safety/quality\\_assurance/guidelines/en/](http://www.who.int/medicines/areas/quality_safety/quality_assurance/guidelines/en/)) for comments under the “Current projects” link. If you wish to receive our draft guidelines, please send your e-mail address to [jonessi@who.int](mailto:jonessi@who.int) and your name will be added to our electronic mailing list.

---

© World Health Organization 2020

All rights reserved.

This draft is intended for a restricted audience only, i.e. the individuals and organizations having received this draft. The draft may not be reviewed, abstracted, quoted, reproduced, transmitted, distributed, translated or adapted, in part or in whole, in any form or by any means outside these individuals and organizations (including the organizations' concerned staff and member organizations) without the permission of the World Health Organization. The draft should not be displayed on any website.

Please send any request for permission to:

Dr Sabine Kopp, Team Lead, Norms and Standards for Pharmaceuticals, Technical Standards and Specifications, Health Products Policy and Standards Department, Access to Medicines and Health Products Division, World Health Organization, CH-1211 Geneva 27, Switzerland, email: [kopps@who.int](mailto:kopps@who.int).

The designations employed and the presentation of the material in this draft do not imply the expression of any opinion whatsoever on the part of the World Health Organization concerning the legal status of any country, territory, city or area or of its authorities, or concerning the delimitation of its frontiers or boundaries. Dotted lines on maps represent approximate border lines for which there may not yet be full agreement.

The mention of specific companies or of certain manufacturers' products does not imply that they are endorsed or recommended by the World Health Organization in preference to others of a similar nature that are not mentioned. Errors and omissions excepted, the names of proprietary products are distinguished by initial capital letters.

All reasonable precautions have been taken by the World Health Organization to verify the information contained in this draft.

However, the printed material is being distributed without warranty of any kind, either expressed or implied. The responsibility for the interpretation and use of the material lies with the reader. In no event shall the World Health Organization be liable for damages arising from its use.

This draft does not necessarily represent the decisions or the stated policy of the World Health Organization.

SCHEDULE FOR THE PROPOSED ADOPTION PROCESS OF DOCUMENT QAS/20.837:  
**DISSOLUTION TEST FOR SOLID ORAL DOSAGE FORMS**

Description	Date
Revision drafted.	February 2020
Discussion at the consultation on Screening Technologies, Laboratory Tools and Pharmacopoeial Specifications for Medicines.	May 2020
Draft revision sent out for public consultation.	July - August 2020
Presentation to the 55 <sup>th</sup> WHO Expert Committee on Specifications for Pharmaceutical Preparations.	October 2020
Further follow-up action as required.	

*[Note from the Secretariat. Chapter 5.5 Dissolution test for solid oral dosage form is based on the internationally-harmonized texts developed by the Pharmacopoeial Discussion Group (PDG). It is proposed to add to the chapter a section for the analysis of Suspensions and powders for suspension, which did not undergo harmonization and is not part of the PDG text. Comments are sought on the new section.*

*Changes from the current text are indicated in the text by insert or ~~delete~~.*

## Revision of chapter 5.5 DISSOLUTION TEST FOR SOLID ORAL DOSAGE FORMS

### 5.5 DISSOLUTION TEST FOR ~~SOLID~~ ORAL DOSAGE FORMS

*This text is based on the internationally-harmonized texts developed by the Pharmacopoeial Discussion Group (PDG). It has been developed in line with the style and requirements used in The International Pharmacopoeia. The additional sections on Suspensions and powders for oral suspension and on monographs of The International Pharmacopoeia ~~are~~ *is not part of the PDG text.**


*For further guidance, see also the chapter Dissolution testing of tablets and capsules in the Supplementary Information section.*

This test determines the amount of active ingredient(s) released from an solid-oral dosage form, such as a tablet or a capsule, under controlled conditions using a known volume of dissolution medium within a predetermined length of time.

**Basket apparatus.** The assembly consists of the following: a vessel, which may be covered, made of glass or other inert, transparent material, which should not sorb, react or interfere with the preparation to be tested; a motor; a drive shaft; and a cylindrical basket (stirring element). The vessel is partially immersed in a suitable water-bath of any convenient size or heated by a suitable device such as a heating jacket. The water-bath or heating device permits maintaining the temperature inside the vessel at  $37 \pm 0.5$  °C during the test. No part of the assembly, including the environment in which the assembly is placed, contributes significant motion, agitation or vibration beyond that due to the smoothly rotating stirring element. An apparatus that permits observation of the preparation and stirring element during the test is preferable. The vessel is cylindrical with a hemispherical bottom and a capacity of 1 litre. Its height is 160–210 mm and its inside diameter is 98–106 mm. Its sides are flanged at the top. A fitted cover may be used to retard evaporation. If a cover is used, it provides sufficient openings to allow ready insertion of the thermometer and withdrawal of samples. The shaft is positioned so that its axis is not more than 2 mm at any point from the vertical axis of the vessel and rotates smoothly and without significant wobble that could affect the results. A speed-regulating

device is used that allows the shaft rotation speed to be selected and maintained at a specified rate within  $\pm 4\%$ .


Shaft and basket components of the stirring element are fabricated of stainless steel, type 316 or equivalent, to the specifications shown in Figure 1. A basket having a gold coating of about  $2.5\ \mu\text{m}$  (0.0001 inch) thick may be used. The dosage unit is placed in a dry basket at the beginning of each test. The distance between the inside bottom of the vessel and the bottom of the basket is maintained at  $25 \pm 2\ \text{mm}$  during the test.

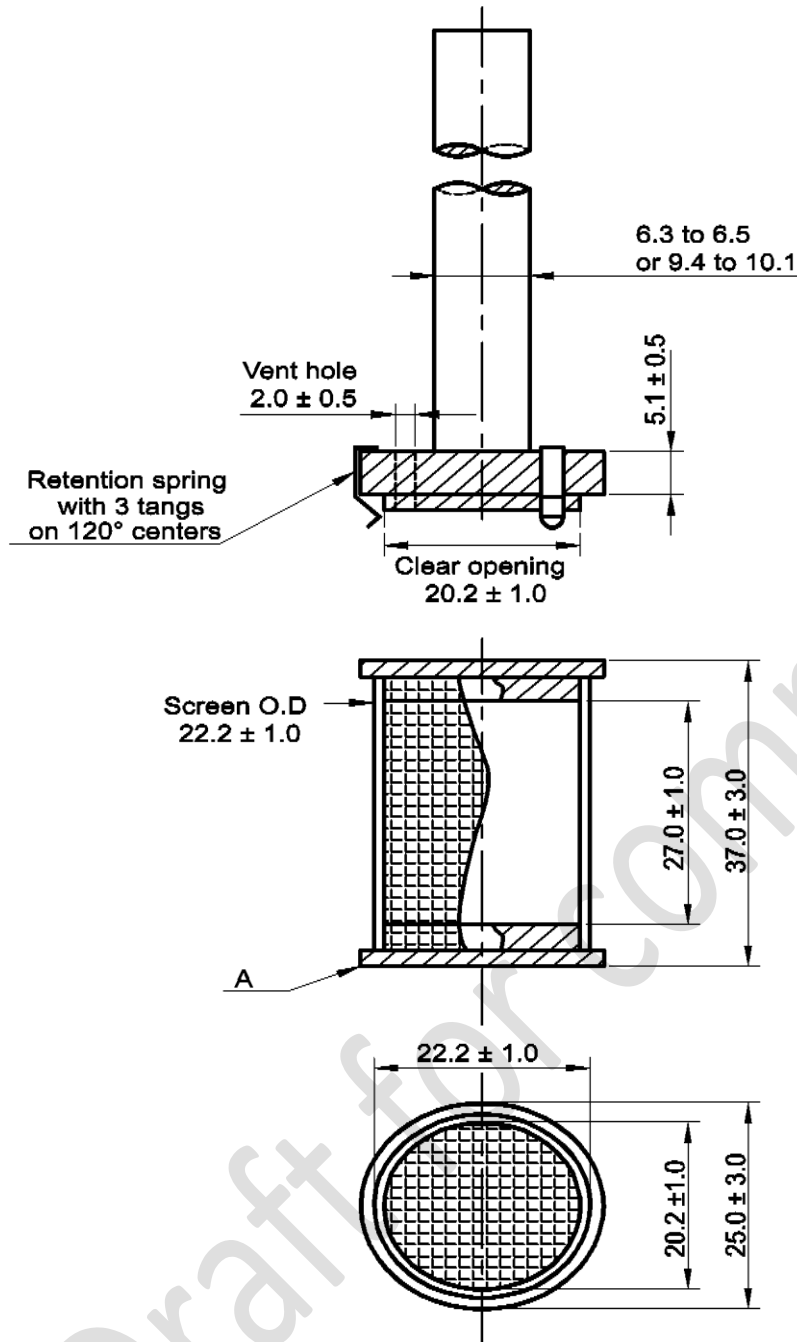

1. Screen with welded seam: 0.22–0.31 mm wire diameter with wire opening of 0.36–0.44 mm. After welding, the screen may be slightly altered.
2. Maximum allowable runout at “A” is 1.0 mm when the part is rotated on centre line axis with basket mounted.

*Figure 1. Basket stirring element*

Dimensions in millimetres.

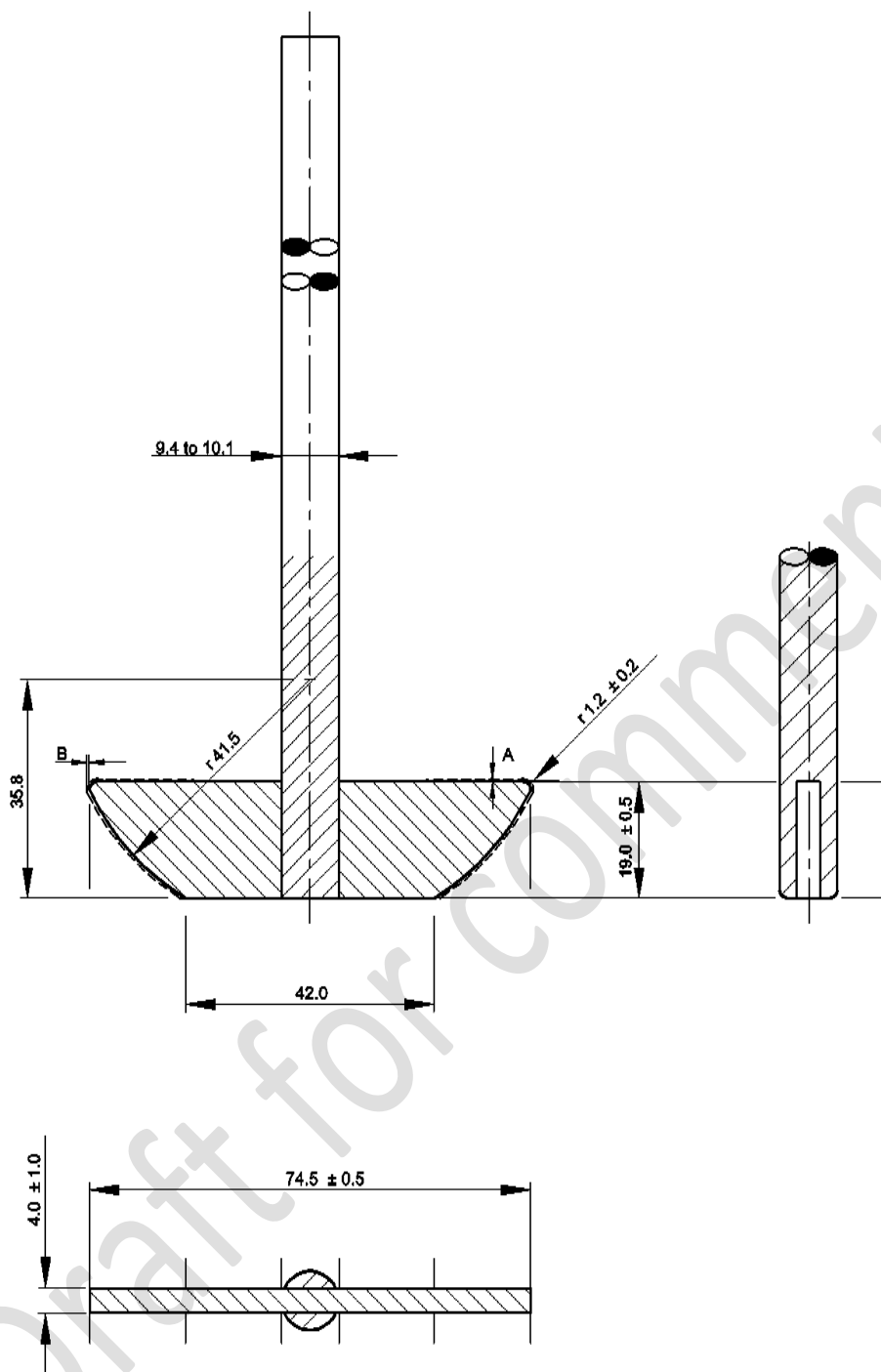

109

110

111 A and B dimensions do not vary more than 0.5 mm when part is rotated on centre line axis.

112 Tolerances are  $\pm 1.0$  mm unless otherwise stated.

113

114 *Figure 2. Paddle stirring element*

115

116 Dimensions in millimetres.

**Paddle apparatus.** Use the assembly from the basket apparatus except that a paddle formed from a blade and a shaft is used as the stirring element. The shaft is positioned so that its axis is not more than 2 mm from the vertical axis of the vessel at any point and rotates smoothly without significant wobble that could affect the results. The vertical centre line of the blade passes through the axis of the shaft so that the bottom of the blade is flush with the bottom of the shaft. The paddle conforms to the specifications shown in Figure 2. The distance of  $25 \pm 2$  mm between the bottom of the blade and the inside bottom of the vessel is maintained during the test. The metallic or suitably inert, rigid blade and shaft comprise a single entity. A suitable two-part detachable design may be used provided the assembly remains firmly engaged during the test. The paddle blade and shaft may be coated with a suitable coating so as to make them inert. The dosage unit is allowed to sink to the bottom of the vessel before rotation of the blade is started. A small, loose piece of non-reactive material, such as not more than a few turns of wire helix, may be attached to dosage units that would otherwise float. An alternative sinker device is shown in Figure 3. Other validated sinker devices may be used.

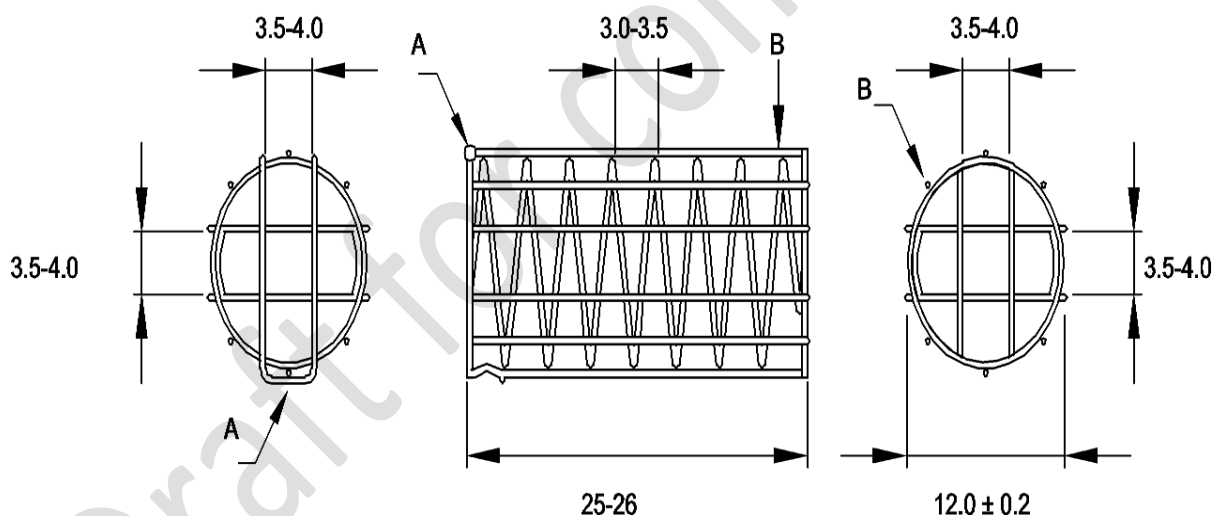


Figure 3. Alternative sinker. A: acid-resistant wire clasp; B: acid-resistant wire support; Dimensions in millimeters.

## Recommended procedure

### *Conventional-release (or immediate-release) dosage forms*

*Procedure.* Place the stated volume of the dissolution medium ( $\pm 1\%$ ) in the vessel of the specified apparatus. Assemble the apparatus, equilibrate the dissolution medium to  $37 \pm 0.5$  °C and remove the thermometer. The test may also be carried out with the thermometer in place, provided it is shown that results equivalent to those obtained without the thermometer are obtained. Place one dosage unit in the apparatus taking care to exclude air bubbles from the surface of the dosage unit. Operate the apparatus at the specified rate. Within the time interval specified, or at each of the times stated, withdraw a sample from a zone midway between the surface of the dissolution medium and the top of the rotating basket or blade not less than 1 cm from the vessel wall. Agitation/stirring should continue during sampling. Where multiple sampling times are specified replace the samples withdrawn for analysis with equal volumes of fresh dissolution medium at 37 °C or, where it can be shown that replacement of the medium is not necessary, correct for the volume change in the calculation. Keep the vessel covered for the duration of the test and verify the temperature ( $37 \pm 0.5$  °C) of the medium at suitable times. Perform the analysis as directed in the individual monograph using a suitable assay method. The samples are filtered immediately upon sampling, preferably by using in-line filtration or a filter in the tip of the sampling probe or both, unless filtration is demonstrated to be unnecessary. Use an inert filter that does not cause adsorption of the active ingredient or contain extractable substances that would interfere with the analysis. Centrifugation is not recommended unless validated for the specific test. The test is to be conducted with six dosage form units in parallel.

If automated equipment is used for sampling or the apparatus is otherwise modified verification is necessary that the modified apparatus will produce results equivalent to those obtained with the apparatus described in this chapter.

*Dissolution medium.* A suitable dissolution medium is used. The volume specified refers to measurements made between 20 °C and 25 °C. If the dissolution medium is a buffered solution adjust the solution so that its pH is within 0.05 units of the specified pH. Dissolved gases can



cause bubbles to form which may change the results of the test. In such cases, dissolved gases must be removed prior to testing.<sup>1</sup>

*Time.* Where a single time specification is given, the test may be concluded in a shorter period if the requirement for minimum amount dissolved is met. Samples are to be withdrawn only at the stated times, within a tolerance of  $\pm 2\%$ .

Determine the quantity of active ingredient dissolved at the specified time(s) indicated in the individual monograph. The result should be expressed as a percentage of the content stated on the label of the dosage form.

#### **Sustained-release (or extended-/prolonged-release) solid-dosage forms**

*Procedure.* Proceed as described for conventional-release dosage forms.

*Dissolution medium.* Proceed as described for conventional-release dosage forms.

*Time.* The test-time points, generally three, are expressed in hours.

#### **Delayed-release solid-dosage forms**

*Procedure.* Use method A or B.

##### **Method A**

- *Acid stage.* Place 750 mL *hydrochloric acid (0.1 mol/L) VS* in the vessel and assemble the apparatus. Allow the medium to equilibrate to a temperature of  $37 \pm 0.5$  °C. Place one dosage unit in the apparatus, cover the vessel and operate the apparatus at the specified rate. After 2 hours of operation in *hydrochloric acid (0.1 mol/L) VS*, withdraw

---

<sup>1</sup> One appropriate method of deaeration is as follows: heat the medium, while stirring gently, to about 41 °C, immediately filter under vacuum using a filter having a pore size of 0.45 µm or less, with vigorous stirring and continue stirring under vacuum for at least 5 minutes, preferably 15 minutes, until no more bubbles are observed. Other validated deaeration techniques for removal of dissolved gases may be used.

a sample of the fluid and proceed immediately as directed under buffer stage. Perform an analysis of the sample using a suitable assay method.

- *Buffer stage.* Complete the operations of adding and adjusting the pH within 5 minutes. With the apparatus operating at the rate specified, add to the fluid in the vessel 250 mL of a 0.2 M solution of *trisodium orthophosphate R* that has been equilibrated to  $37 \pm 0.5$  °C. Adjust, if necessary, with *hydrochloric acid (~70 g) TS* or *sodium hydroxide (~80 g/L) TS* to a pH of  $6.8 \pm 0.05$ . Continue to operate the apparatus for 45 minutes or for the specified time. At the end of the time period, withdraw a sample of the fluid and perform the analysis using a suitable assay method.

#### Method B

- *Acid Stage.* Place 1000 mL of *hydrochloric acid (0.1 mol/L) VS* in the vessel and assemble the apparatus. Allow the medium to equilibrate to a temperature of  $37 \pm 0.5$  °C. Place one dosage unit in the apparatus, cover the vessel and operate the apparatus at the specified rate. After 2 hours of operation in *hydrochloric acid (0.1 mol/L) VS*, withdraw a sample of the fluid and proceed immediately as directed under buffer stage. Perform an analysis of the sample using a suitable assay method.

- *Buffer stage.* For this stage of the procedure, use buffer that has previously been equilibrated to a temperature of  $37 \pm 0.5$  °C. Drain the acid from the vessel and add 1000 mL of pH 6.8 phosphate buffer, prepared by mixing three volumes of *hydrochloric acid (0.1 mol/L) VS* with one volume of a 0.20 M solution of *trisodium orthophosphate R* and adjusting, if necessary, with *hydrochloric acid (~70 g/L) TS* or *sodium hydroxide (~80 g/L) TS* to a pH of  $6.8 \pm 0.05$ . This may also be accomplished by removing from the apparatus the vessel containing the acid and replacing it with another vessel containing the buffer and transferring the dosage unit to the vessel containing the buffer. Continue to operate the apparatus for 45 minutes or for the specified time. At the end of the time period, withdraw a sample of the fluid and perform the analysis using a suitable assay method.

*Time.* All test times stated are to be observed within a tolerance of  $\pm 2\%$ , unless otherwise specified.

## Acceptance criteria

### *Conventional-release (or immediate-release) dosage forms*

Unless otherwise specified in the individual monograph, the requirements are met if the quantities of active ingredient(s) dissolved from the dosage forms tested conform to Table 1. Continue testing through the three levels unless the results conform at either  $S_1$  or  $S_2$ . The quantity,  $Q$ , is the specified amount of dissolved active ingredient expressed as a percentage of the labelled content; the 5%, 15% and 25% values in the acceptance table are percentages of the labelled content so that these values and  $Q$  are in the same terms.

Table 1

Level	Samples tested	Acceptance criteria
$S_1$	6	Each value is not less than $Q + 5\%$ .
$S_2$	6	Average value of the 12 dosage units ( $S_1 + S_2$ ) is equal to or greater than $Q$ and no unit is less than $Q - 15\%$ .
$S_3$	12	Average value of 24 dosage units ( $S_1 + S_2 + S_3$ ) is equal to or greater than $Q$ ; not more than 2 units are less than $Q - 15\%$ ; no unit is less than $Q - 25\%$ .

### **Sustained-release (or extended-/prolonged-release) dosage forms**

Unless otherwise specified in the individual monograph, the requirements are met if the quantities of active ingredient(s) dissolved from the dosage forms tested conform to Table 2. Continue testing through the three levels unless the results conform at either  $L_1$  or  $L_2$ . Limits on the amounts of active ingredient(s) dissolved are expressed in terms of the labelled content. The limits embrace each value of  $Q_i$ , the amount dissolved at each specified fractional dosing interval. Where more than one range is specified, the acceptance criteria apply individually to each range.

Table 2

Level	Samples tested	Acceptance criteria
L <sub>1</sub>	6	No individual value lies outside each of the stated ranges and no individual value is less than the stated amount at the final test time.
L <sub>2</sub>	6	The average value of the 12 dosage units (L <sub>1</sub> + L <sub>2</sub> ) lies within each of the stated ranges and is not less than the stated amount at the final test time; none is more than 10% of the labelled content outside each of the stated ranges; and none is more than 10% of labelled content below the stated amount at the final test time.
L <sub>3</sub>	12	The average value of the 24 dosage units (L <sub>1</sub> + L <sub>2</sub> + L <sub>3</sub> ) lies within the stated ranges and is not less than the stated amount at the final test time; not more than 2 of the 24 dosage units are more than 10% of labelled content outside each of the stated ranges; not more than 2 of the 24 dosage units are more than 10% of labelled content below the stated amount at the final test time; and none is more than 20% of labelled content outside each of the stated ranges or more than 20% of labelled content below the stated amount at the final test time.

### Delayed-release dosage forms

*Acid stage.* Unless otherwise stated in the individual monograph, the requirements of this part of the test are met if the quantities, based on the percentage of the labelled content of active ingredient(s) dissolved from the dosage units tested conform to Table 3. Continue testing through the three levels unless the results of both acid and buffer stages conform at an earlier level.

Table 3

Level	Samples tested	Acceptance criteria
A <sub>1</sub>	6	No individual value exceeds 10% dissolved.
A <sub>2</sub>	6	Average value of the 12 dosage units (A <sub>1</sub> + A <sub>2</sub> ) is not more than 10% dissolved, and no individual value is greater than 25% dissolved.

A <sub>3</sub>	12	Average value of 24 dosage units ( $A_1 + A_2 + A_3$ ) is not more than 10% dissolved, and no individual value is greater than 25% dissolved.
----------------	----	---

*Buffer stage.* Unless otherwise specified in the individual monograph, the requirements are met if the quantities of active ingredients dissolved from the units tested conform to Table 4. Continue testing through the three levels unless the results of both stages conform at an earlier level. The value of  $Q$  in Table 4 is 75% dissolved unless otherwise specified. The quantity,  $Q$ , is the specified total amount of active ingredient dissolved in both the acid and buffer stages, expressed as a percentage of the labelled content. The 5%, 15% and 25% values in the table are percentages of the labelled content so that these values and  $Q$  are in the same terms.

Table 4

Level	Samples tested	Acceptance criteria
B <sub>1</sub>	6	No value is less than $Q + 5\%$ .
B <sub>2</sub>	6	Average value of the 12 dosage units ( $B_1 + B_2$ ) is equal to or greater than $Q$ , and no unit is less than $Q - 15\%$ .
B <sub>3</sub>	12	Average value of the 24 dosage units ( $B_1 + B_2 + B_3$ ) is equal to or greater than $Q$ ; not more than 2 units are less than $Q - 15\%$ , and no unit is less than $Q - 25\%$ .

### **Suspensions and Powders for oral suspension**

*Procedure.* Prepare a suspension from Powders for oral suspension or re-suspend suspensions according to the product information. Avoid the introducing of bubbles into the sample to ensure the precision of dosing. Transfer to the dissolution apparatus, preferably by weight, the amount of sample indicated in the monograph. If no amount is given, use an amount that is equivalent to 1 unit dose or, in case the product has different doses depending on body weight or age, the amount of sample that corresponds to the highest unit dose to be administered at one time. If the product is labeled for single use, each sample should come from a different container/packet.

## **Monographs of *The International Pharmacopoeia***

*The following additional statements apply to the individual monographs of The International Pharmacopoeia.*

### **Qualification of dissolution test equipment and verification of system performance<sup>2</sup>**

Periodically qualify the equipment utilizing an “enhanced mechanical calibration”, such as the procedure described in the international standard procedure ASTM 2503 or a combination of a mechanical calibration to determine conformance of the dissolution apparatus to the dimensions and tolerances as given above and the analysis of suitable reference tablets to verify the performance of the testing system.

### **Test conditions**

The following specifications are given in the individual monographs:

- the apparatus to be used;
- the composition and volume of the dissolution medium;
- the rotation speed of the paddle or basket;
- the preparation of the test and reference solutions;
- the time, the method and the amount of sample to be withdrawn or the conditions for continuous monitoring; the preparation of the sample and the reference solution;
- the method of analysis; and
- the limits of the quantity or quantities of active pharmaceutical ingredient(s) required to dissolve within a prescribed time.

### **Dissolution media**

If a buffer is added to the dissolution medium, adjust its pH to within  $\pm 0.05$  units of the prescribed value.

---

<sup>2</sup> See also Supplementary Guidelines on Good Manufacturing Practices: Validation. World Health Organization. WHO Technical Report Series, No. 937, 2006.

In specific cases, and subject to approval by the relevant regional or national authority, dissolution media may contain enzymes and/or surfactants. The addition of enzymes may be considered, for example, for formulations containing gelatin when dissolution failures can be ascribed to the cross-linking of this excipient (e.g. hard and soft gelatin capsules; gelatin containing tablets). For the testing of preparations containing poorly aqueous-soluble active substances, modification of the medium may be necessary. A surfactant may be added only when the active pharmaceutical ingredient is insoluble over the entire physiological pH range, pH 1.2 to 6.8. In such circumstances, a low concentration of surfactant may be prescribed.

Below are some examples of dissolution media:

- **Dissolution buffer pH 1.2, TS**

Dissolve 2 g of *sodium chloride R* in 800 mL of *water R*, adjust the pH to 1.2 with *hydrochloric acid (~70 g/L) TS* and dilute to 1000 mL with *water R*.

- **Dissolution buffer pH 2.5, TS**

Dissolve 2 g of *sodium chloride R* in 800 mL of *water R*, adjust the pH to 2.5 with *hydrochloric acid (~70 g/L) TS* and dilute to 1000 mL with *water R*.

- **Dissolution buffer pH 3.5, TS**

Dissolve 7.507 g of *glycine R* and 5.844 g of *sodium chloride R* in 800 mL of *water R*, adjust the pH to 3.5 with *hydrochloric acid (~70 g/L) TS* and dilute to 1000 mL with *water R*.

- **Dissolution buffer pH 4.5, TS1**

Dissolve 2.99 g of *sodium acetate R* in 900 mL of *water R*, adjust the pH to 4.5 by adding about 14 mL of *acetic acid (~120 g/L) TS* and dilute to 1000 mL with *water R*.

• **Dissolution buffer pH 4.5, TS2**

Dissolve 6.8 g of *potassium dihydrogen phosphate R* in 900 mL of *water R*, adjust the pH to 4.5 either with *hydrochloric acid (~70 g/L) TS* or *sodium hydroxide (~80 g/L) TS* and dilute to 1000 mL with *water R*.

• **Dissolution buffer, pH 6.8, TS**

Dissolve 6.9 g of *sodium dihydrogen phosphate R* and 0.9 g of *sodium hydroxide R* in 800 mL of *water R*, adjust the pH to 6.8 with *sodium hydroxide (~80g/L) TS* and dilute to 1000 mL with *water R*.

• **Dissolution buffer, pH 6.8, 0.25% SDS TS**

Dissolve 6.9 g of *sodium dihydrogen phosphate R*, 0.9 g of *sodium hydroxide R* and 2.5 g of *sodium dodecyl sulfate R* in 800 mL of *water R*, adjust the pH to 6.8 with *sodium hydroxide (~80g/L) TS* and dilute to 1000 mL with *water R*.

• **Dissolution buffer pH 7.2, TS**

Dissolve 9.075 g of *potassium dihydrogen phosphate R* in *water R* to produce 1000 mL (solution A). Dissolve 11.87 g of *disodium hydrogen phosphate R* in sufficient *water R* to produce 1000 mL (solution B). Mix 300 mL of solution A with 700 mL of solution B.

• **Gastric fluid, simulated, TS**

Dissolve 2.0 g of *sodium chloride R* and 3.2 g of *pepsin R* in 7.0 mL of *hydrochloric acid (~420 g/L) TS* and sufficient *water R* to produce 1000 mL. This test solution has a pH of about 1.2.



• **Intestinal fluid pH 6.8, simulated, TS**

Mix 77.0 mL of *sodium hydroxide* (0.2 mol/L) VS, 250.0 mL of a solution containing 6.8 g *potassium dihydrogen phosphate R* and 500 mL of *water R*. Add 10.0 g *pancreatin R*, mix and adjust the pH with the buffer components to  $6.8 \pm 0.1$ . Dilute to 1000 mL with *water R*.

\*\*\*

Draft for comments